



JPES Journal of Physical Education and Sport



Online Publication Date: 10 December 2009

ORIGINAL RESEARCH

SURFACE ELECTROMYOGRAPHY IN BIOMECHANICS: APPLICATIONS AND SIGNAL ANALYSIS ASPECTS

DEAK GRAȚIELA-FLAVIA¹, **RUSU FLAVIA¹**, GROSU EMILIA¹¹Babeș-Bolyai University of Cluj-Napoca

Physical Education and Sport Faculty, gratiela_flavia@yahoo.com

Abstract. Surface electromyography (SEMG) is a technique for detecting and recording the electrical activity of the muscles using surface electrodes. The EMG signal is used in biomechanics mainly as an indicator of the initiation of muscle activation, as an indicator of the force produced by a contracting muscle, and as an index of the fatigue occurring within a muscle. EMG, used as a method of investigation, can tell us if the muscle is active or not, if the muscle is more or less active, when it is on or off, how much active is it, and finally, if it fatigues. The purpose of this article is to discuss some specific EMG signal analysis aspects with emphasis on comparison type analysis and frequency fatigue analysis.

Keywords: electromyography, signal, biomechanics, applications, sport

REZUMAT. *Electromiografia de suprafață în biomecanică: aplicații și aspecte ale analizei de semnal.* Electromiografia de suprafață (SEMG) este o tehnică de detecție și înregistrare a activității electrice a mușchilor cu ajutorul electrozilor de suprafață. Semnalul electromiografic este folosit în biomecanică în principal ca un indicator al inițierii activării musculare, ca un indicator al forței produse de un mușchi în contracție și, nu în ultimul rând, ca un index al oboselii musculare. Electromiografia, folosită ca metodă de investigație, ne poate spune dacă mușchiul este activ sau nu, dacă este mai mult sau mai puțin activ, când este activ și când nu, cât de mult este activ și dacă este obosit. Scopul acestui articol este acela de a discuta câteva aspecte specifice analizei de semnal electromiografic, punând accentul pe analizele comparative și pe cele de frecvență-oboseală.

CUVINTE CHEIE: electromiografie, semnal, biomecanică, aplicații, sport

Introduction

Electromyography (EMG) is a technique for detecting and recording the electrical activity of the muscles (Strungaru, 1982). The instrument used to detect the electrical potential generated by the muscle cells when they contract or are at rest is called an electromyograph. The record produced by an electromyograph is known as an electromyogram or a myogram. The signal recorded as a myogram is the result of superimposed motor unit action potentials (MUAPs) from several motor units. A motor unit has, as a central element, one motor neuron which innervates several muscle fibers. The impulse generated by a firing motor unit is called an action potential and it is carried down the motor neuron to the muscle. The neuromuscular junction is the area

where the nerve contacts the muscle. Once the action potential is transmitted across the neuromuscular junction, an action potential is elicited in all of the innervated muscle fibers of that particular motor unit. The sum of all this electrical activity is known as a *motor unit action potential* (MUAP) (Wikipedia).

The raw EMG signal looks like the one presented in Figure 1. The amplitude of the signal can range from 0 to 10 mV (peak-to-peak). The usable energy of the signal ranges between 0 and 500 Hz, with the dominant energy being in the 50-150 Hz range (De Luca, 2002). When analyzed, the EMG signal is rarely used as a raw one. Typically, the information derived from a raw EMG signal is limited to a frequency analysis. In order to interpret the amplitude and time characteristics, the EMG signal must be processed. In most cases, the steps followed in processing the raw EMG signal are rectification, smoothing and digital filtering (Konrad, 2005).

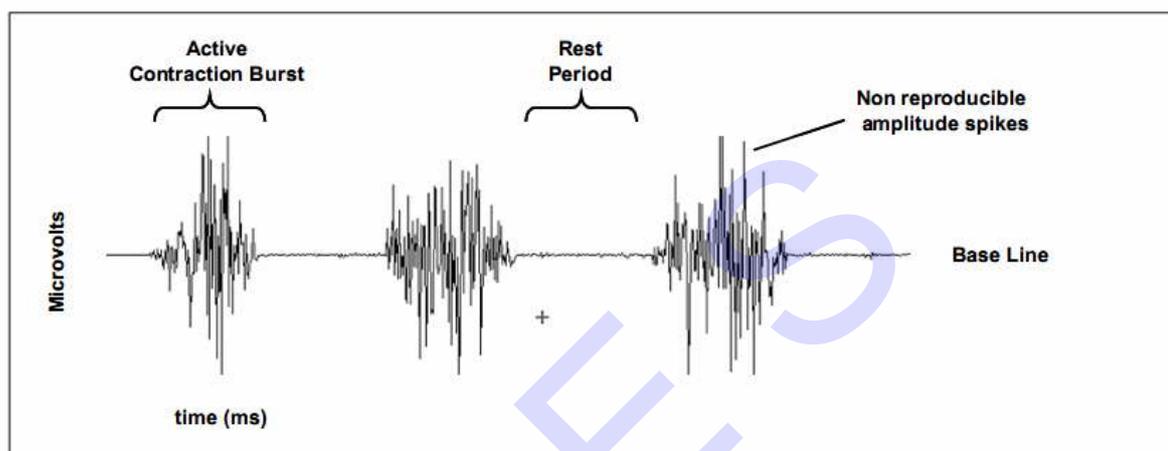


Figure 1. The raw EMG signal (Adopted from Konrad, 2005)

Surface electromyography (SEMG) is mainly used in kinesiological studies because of its non-invasive character. It implies the detection of the myoelectric signal with surface electrodes. The main drawback of SEMG is the fact that only surface muscles can be investigated.

Electromyography was initially used in neurology as an analyzing method of the muscle response at an artificial electrical stimulus. Nowadays, EMG is used in medical research, ergonomics, rehabilitation and sports science. In kinesiology, EMG is described as “the study of the neuromuscular activation of muscles within postural tasks, functional movements, work conditions and treatment/training regimes” (Konrad, 2005).

SEMG in biomechanics

In the world of biomechanics, electromyography has the role of evaluating the neuromuscular activation within any type of physical activity. Aside Anthropometry, Kinematics and Kinetics, Kinesiological EMG (Figure 2) is recognized as an objective measurement method (Konrad, 2005). The aims of kinesiological EMG are to analyze the function and coordination of muscles during any kind of physical performance, in both healthy and disabled subjects.

The research areas of kinesiological EMG include studies of normal muscle function during different types of contraction; studies of muscle activity in sports, occupational activity and rehabilitation; evaluation of functional anatomical muscle activity; coordination and synchronization studies; studies on how training affects muscle function; fatigue studies; muscle force and EMG; the influence of equipment on muscle activity, and so on (Clarys, 2000).

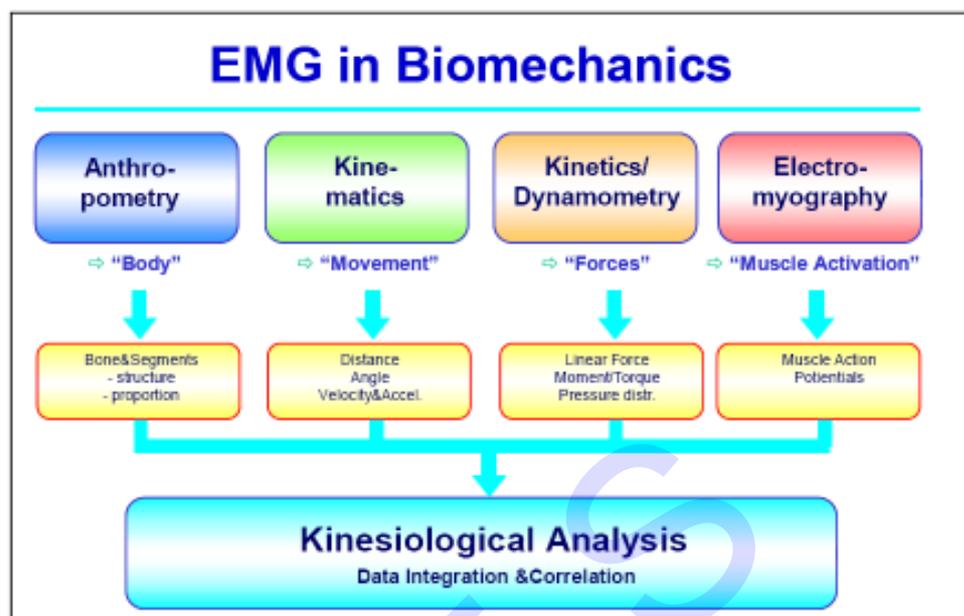


Figure 2. The 4 major areas of biomechanical measurement methods (Adopted from Konrad, 2005)

Surface EMG signal is used in biomechanics mainly as an indicator of the initiation of muscle activation, as an indicator of the force produced by a contracting muscle, and as an index of the fatigue occurring within a muscle. When used as an indicator of the initiation of muscle activity, the SEMG signal provides the timing sequence of one or more muscles which contract with the purpose of fulfilling a task, such as bending forward from an upright position or maintaining the upright standing position (De Luca, 1997). Although the amplitude of the EMG signal does not necessarily reflect the value of the force generated by a contracting muscle, it is qualitatively related to the amount of torque (or force) measured about a joint. This is the reason why some amount of information regarding muscle force can be extracted from EMG data, but this process must be approached with caution (Royer, 2005). If we are interested in assessing the muscle fatigue, the EMG signal is a good information source. The frequency spectrum of the EMG signal, when correctly interpreted, reveals when the investigated muscles fatigue.

SEMG analysis issues

As we have seen so far, EMG focuses on muscles. To be more specific, it focuses on the electrical activity of the muscles. What exactly can EMG, used as a method of investigation, tell us about muscles? First and most important, EMG can tell us if the muscle is active or not. Then it can tell us if the muscle is more or less active, when it is on or off, how much active is it, and finally, if it fatigues (Konrad, 2005). In the following pages we will focus only on three of these issues. We will discuss why it is relevant to know if the muscle is active or not, if the muscle is more or less active, and if it fatigues.

The first question that EMG addresses ("Is the muscle active?") is a simple one. It takes no more than one minute for an EMG specialist to record the raw signal generated by a specific muscle and, thus, respond with Yes or No to the question. A more complex question would be "Is the muscle more or less active?". The answer implies comparison conditions, so the steps to be followed are preparation for recording, recording the signals, analysis of each recorded signal separately, and comparison of the analyzed recordings. Maybe the most complex question that EMG addresses is "Does the muscle fatigue?". As opposite findings were associated with fatigue tests, interpreting the fatigue analysis results is not an easy task. Fulfilling this task requires a great amount of expertise in the field.

There are five major categories of analysis question that EMG can answer. The first one, as we have already mentioned, is referring to whether the muscle is active or not. The answers are yes/no or on/off and they can be obtained by simply observing the raw EMG signal generated by the muscles. The flexion-relaxation phenomenon (Figure 3) exemplifies perfectly the on/off succession of the back muscles when bending forward

from an upright position (Konrad, 2005). At the most flexed position of the trunk, the back muscles (ch. 1 and 2 lumbar right and left; ch. 3 and 4 multifidii right and left) turn off. When slowly extending back, both muscles start firing again.

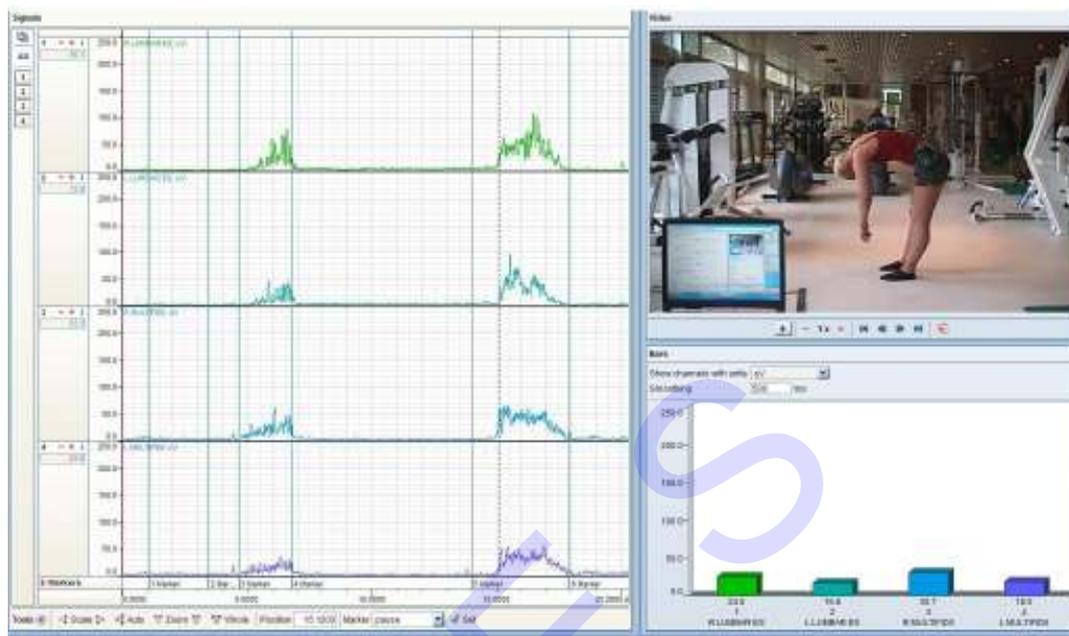


Figure 3. The flexion-relaxation phenomenon (Noraxon MyoResearch XP)

Another example could be the one presented in Figure 4, where a male subject performs 3 bicep curls with his right arm. The signal from the left biceps brachii is recorded on channel 1, and the signal from the right biceps brachii is recorded on channel 2. It can be clearly seen that no signal is recorded on channel 1 as the left biceps brachii is relaxed.

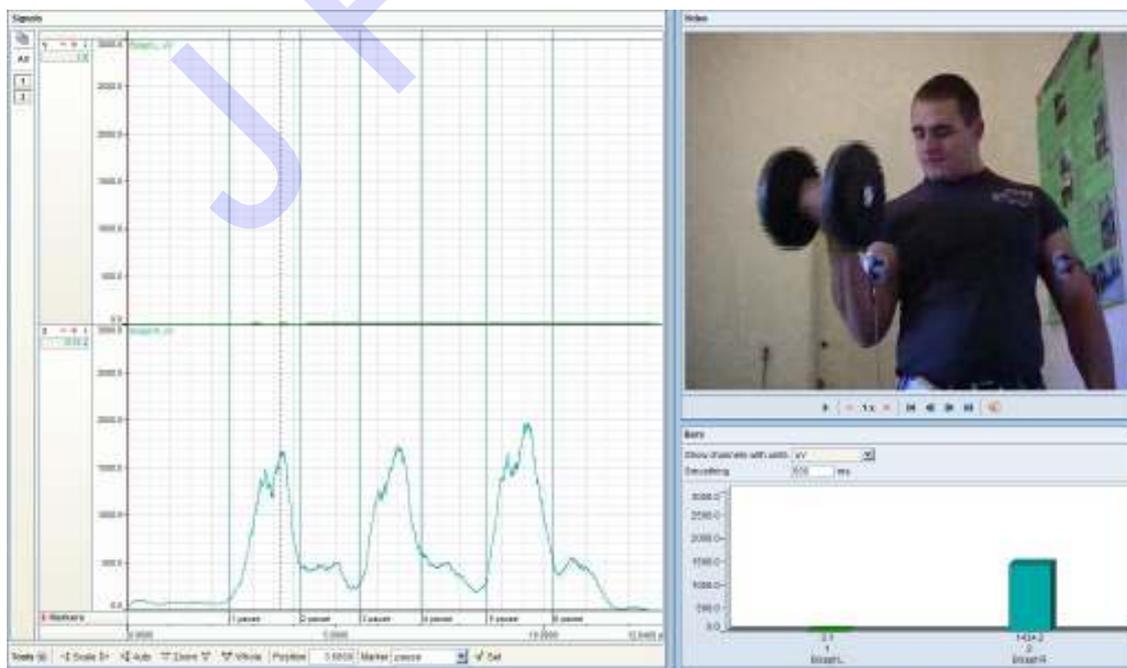


Figure 4. Biceps brachii flexion/extension (Noraxon MyoResearch XP)

At this point, one can argue: we know if one specific muscle is active or not at a given moment. What is the scientific relevance of that? A clinician can and will give a diagnosis based on the results of such a simple EMG analysis. A physical therapist can be provided with a starting point for a therapy regime by simply measuring which muscle of the patient is active in a given task. A training professional can make pertinent choices on the training methods he/she will use to achieve a specific goal by knowing the on/off succession of the muscles in an exercise. These are a few reasons why knowing if a muscle is active or not could be relevant.

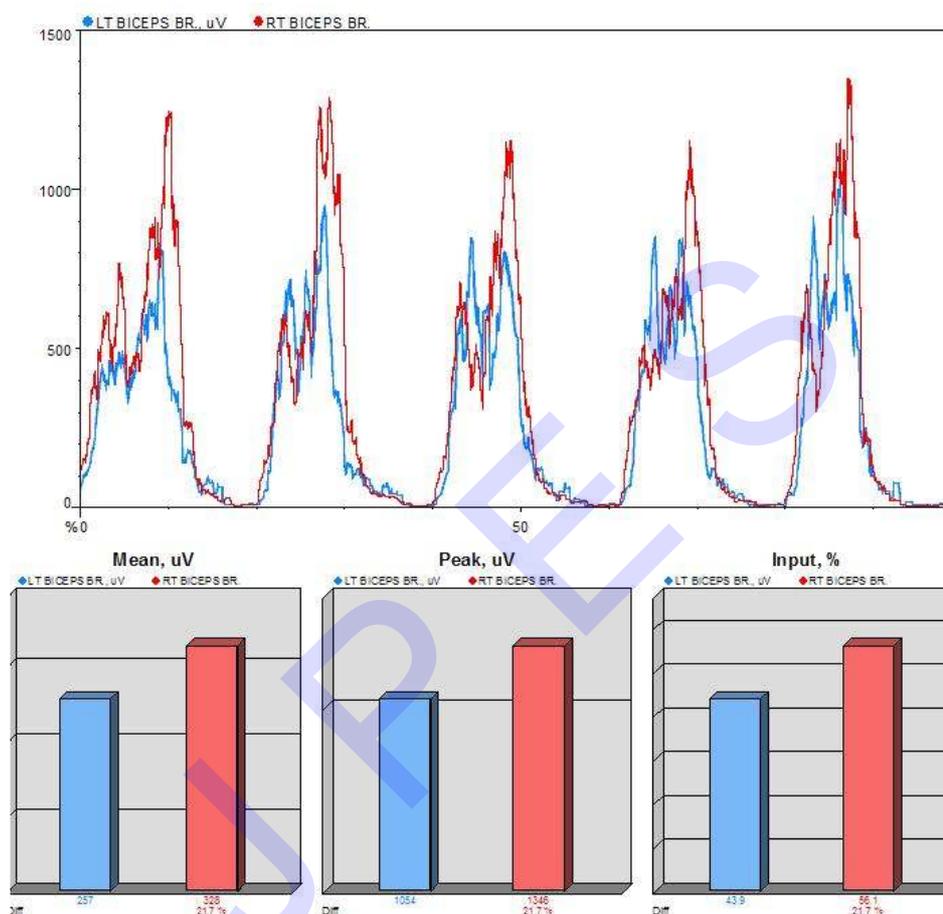


Figure 5. EMG analysis of the biceps brachii in left/right comparison (Noraxon MyoResearch XP)

The question “Is the muscle more or less active?” implies comparison conditions, like comparing the left and right side electromyogram, the pre and post test electromyogram, the subject’s and the normative electromyogram. The answer at this kind of question will be a qualitative one and the data resulted from the qualitative analysis will be presented as curves (Konrad, 2005). A more/less analysis is presented in Figure 5, where the curve with the higher amplitude values is the signal generated by the right biceps brachii of a subject performing 5 voluntary biceps curls, and the curve with the lower amplitude values corresponds to the left biceps brachii of the same subject. Symmetry reports are designed for records measuring the muscle activity in symmetrical, bilateral movements. Comparisons between each muscle pair (left vs. right) are based on EMG signals that are automatically rectified and smoothed online. The MyoResearch XP software’s report overlays the EMG patterns of each muscle pair for visual comparison of pattern, timing and amplitude. Corresponding statistics (mean and peak amplitudes) including the relative difference between each muscle pair are provided in both a bar graph and numerical format (MyoResearch XP, 2004). For the case analyzed in Figure 5, the right biceps brachii is clearly more active than the left one while performing the same task.

Another type of test to which the “more/less” analysis can address is the one presented in Figure 6. This is an ergonomics test and it takes place in an working environment. The muscles tested while the subject performs specific working activities are the lumbar region of the erector spinae (ch. 1), the multifidii (ch. 2), the internal oblique (ch. 3), and the upper trapezius (ch. 4).

